

A man in a dark blue t-shirt and light-colored shorts is lying on his side on a white Ropox Mobilio changing bed. The bed has a grey mattress and white side rails. It is on a light-colored wooden floor against a white wall. A blue text box is overlaid on the left side of the image.

"A movable changing bed gives you greater flexibility"

## ROPOX MOBILIO CHANGING BED

### FREESTANDING AND STABLE

Ropox Mobilio changing bed is mobile and with brake wheels. This ensures that the changing bed can be stable and is easy to move if needed. It is also available with central brake, which locks all wheels at once.

The changing bed can be connected to a standard electrical outlet with a three meter long electrical height adjustment cord. The height is adjustable from 52-92 cm by using a hand control.

A batteri is available as an accessory.

Ropox Mobilio changing bed is available in five models with bed guards on both sides and with a mattress for increased comfort.

The replaceable mattress has a cleaning-friendly surface that makes it easy to clean after use.



**ROPOX™**

BETTER WAYS TO BETTER DAYS






- A PART OF AddLife®

# ROPOX MOBILIO CHANGING BED

Height adjustment: 40 cm (52-92 cm)

Max. load: cf. DS/ISO 17966:2016  
Max. user weight: 200 kg

## VARIANT/MODEL

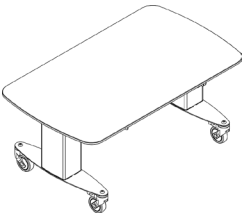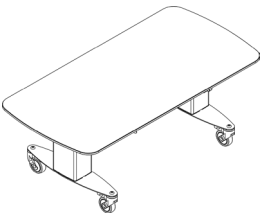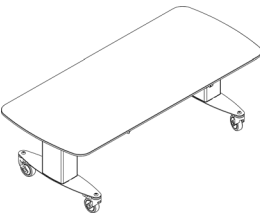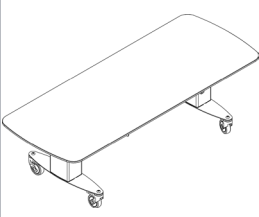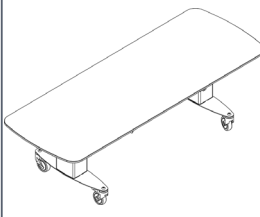MOVABLE CHANGING BED				
				
120x70 cm Item no.: 40-30404	140x70 cm Item no.: 40-30406	160x70 cm Item no.: 40-30408	180x70 cm Item no.: 40-30410	190x70 cm Item no.: 40-30411



## ACCESSORY

Mattress - thickness 3 cm


Length 120 cm Item no.: 40-25760

Length 140 cm Item no.: 40-25795

Length 160 cm Item no.: 40-25796

Length 180 cm Item no.: 40-25764


Length 180 cm Item no.: 40-25765

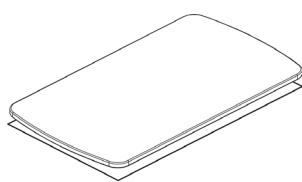

Bed guard

Length 120/140 cm Item no.: 40-30196


Length 160/180/190 cm Item no.: 40-30197

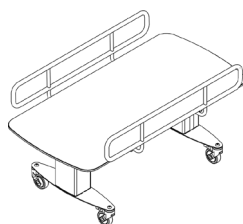

Cover for bed guard

Length 120/140 cm Item no.: 40-25782


Length 160/180/190 cm Item no.: 40-25783

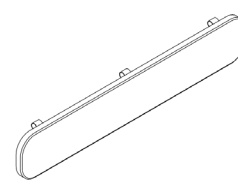

Central brake


Length 120/140 cm Item no.: 40-30765

Length 160/180/190 cm Item no.: 40-30764

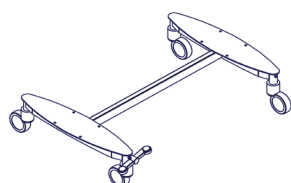

Battery

Item no.: 40-30762

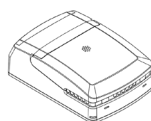

**NB!** Find more information about measurements, technical data, installation interfaces and much more under the product and in the instruction manuals on our website.

The items are produced according to the following Standard: DS/ISO 17966:2016

ROPOX A/S · Ringstedgade 221 · 4700 Naestved, Denmark · Phone: +45 5575 0500 · info@ropox.com · www.ropox.com

**ROPOX™**

BETTER WAYS TO BETTER DAYS

— A PART OF AddLife

No.UK\_1962414\_0123