



## KITCHENS FOR INDEPENDENCE

ACCESSIBILITY AND  
FREEDOM IN THE KITCHEN



## CONTENTS

- 04 THE MULTI-FUNCTIONAL KITCHEN  
We help to choose the right kitchen.
- 06 WORKING AT THE MOST COMFORTABLE HEIGHTS  
Adjust the worktop height to suit the kitchen activity.
- 08 WALL CABINETS THAT COME TO YOU  
We make sure everything is within easy reach.
- 10 ADVICE AND GUIDANCE ON KITCHEN LAYOUTS  
There is more to Ropox than kitchen products.

”Enjoy life in your kitchen – no matter if you are tall or short, sit or stand”

### BECAUSE EVERYDAY LIFE IS IMPORTANT

Cooking and daily kitchen activities are important ingredients of our lives.

**It is our goal that everyone can be free, independent, and safe in their own kitchen.**

For more than 40 years, users around the world have enjoyed the flexibility and accessibility of Ropox kitchen systems. They give the freedom to do most activities independently. With our multi-functional kitchens, more people can cook together or share facilities and be part of the social activities that are so important in our daily lives.

Greater activity improves quality of life.



## The multifunctional kitchen

### ROPOX HELPS TO SELECT THE BEST KITCHEN

We are ready to help and guide you. Together, we can plan a kitchen with optimal space utilization. Involve us from the moment you plan your house and utilize our many years of experience. Together, we match your needs and make sure that the end result is perfect.


Multi-functional kitchens are flexible and offer a wide range of individual solutions for sheltered housing, special schools, nursing homes, training and institutional kitchens, private homes etc. All height-adjustable systems can be hidden and operated either electrically or manually. Easy to fit, the systems are stable and safe and made only of top quality materials.



Ideal both as dining and work table.



Custom-made solutions can be produced.

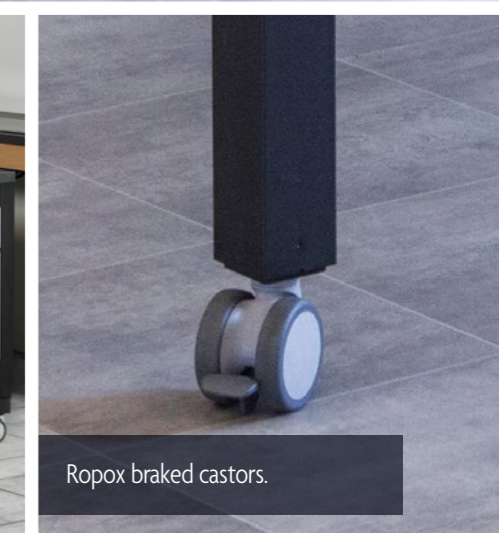

Ropox braked castors.

### ROPOX 4Single frame

- Ideal for the whole family – adults, children and wheelchair users.
- Can be used for a multitude of activities in the kitchen like computer work or art as well as food preparation.
- Three height ranges of 300mm (605-905, 675-975 og 805-11105 mm).
- Adjust the height simply by the push of a button or by using a winding handle or remote control.
- We make the frame to the size of any kitchen worktop or dining tabletop - up to 3000 x 2000mm.
- Hobs and sinks can be integrated as needed.
- Optional braked castors make it easy to move the table, but standard adjustable feet or hidden floor fixings can be used.
- A rechargeable battery solution means the electric version is easily mobile.
- Fully charged batteries can last up to 3 months.
- The frame is designed to withstand a load of up to 150kg (250kg heavy version).



# Ropox worktop systems

## THE WORKTOP REACHES NEW HEIGHTS

Having the worktop adjusted to the perfect height, makes a big difference for each activity. Different heights are needed to knead dough, stir pots, chop vegetables or wash up.

It's a huge advantage to be able to adjust the worktop height – whether standing or seated, promoting self-sufficiency and optimal working conditions in the kitchen. What's more, our solutions offer seated users more leg space and the freedom to move around with ease.




Two separate worktops giving two different heights in the kitchen for two people.

Safety stop to increase safety and to avoid the risk of trapping



## ROPOX FlexiElectric/FlexiManual/FlexiCorner

- Ideal for private homes, schools, training kitchens, sheltered housing, nursing homes etc.
- All types of kitchen worktops can be fitted to the frames in lengths from 600 – 3000mm.
- FlexiCorner suits single person use and can be any size up to 3000x3000mm.
- Sinks and hobs can be integrated in the worktop.
- Easy to adjust in height, manually or electrically.
- Corner solution for both internal and external corners.
- For added safety we always recommend safety stops for electric systems to avoid the risk of trapping.
- RopoxConnect App for adjusting the table.

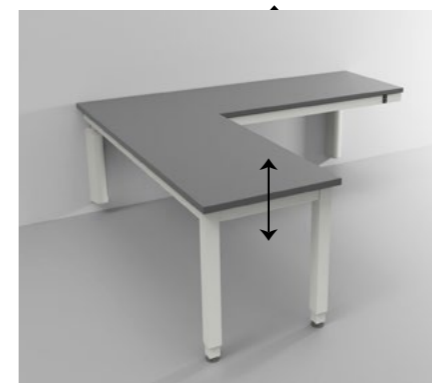

## ROPOX Flexi4Single

- Ideal for private homes, schools, training kitchens, sheltered housing, nursing homes etc.
- Perfect solution for rooms with only one load-bearing wall.
- All types of kitchen worktops can be fitted to the frames.
- Sinks and hobs can be integrated in the worktop.
- Easy to adjust in height, electric with a press pad.
- We always recommend safety stops for electric systems to avoid the risk of trapping.



## ROPOX Worktop brackets

- For places with only rare requirements for height adjustment.
- Requires two people to lift and adjust the height.




# Ropox wall cabinets

## WALL CABINETS THAT COME TO YOU

Our goal is to make it easy to see and reach everything in all wall cabinets – the solution you choose, will depend on ability to reach and if standing or seated. Maybe straight up/down movement brings everything into reach. To make things even easier the Diagonal brings the whole cabinet to the front of the worktop moving over the kettle, toaster or tap. A simple push of a button will adjust the height.

Safety stop plates beneath the units ensure nothing is trapped between the unit and worktop. Again, the needs and situation decide what's best. Together, we can find the perfect solution!

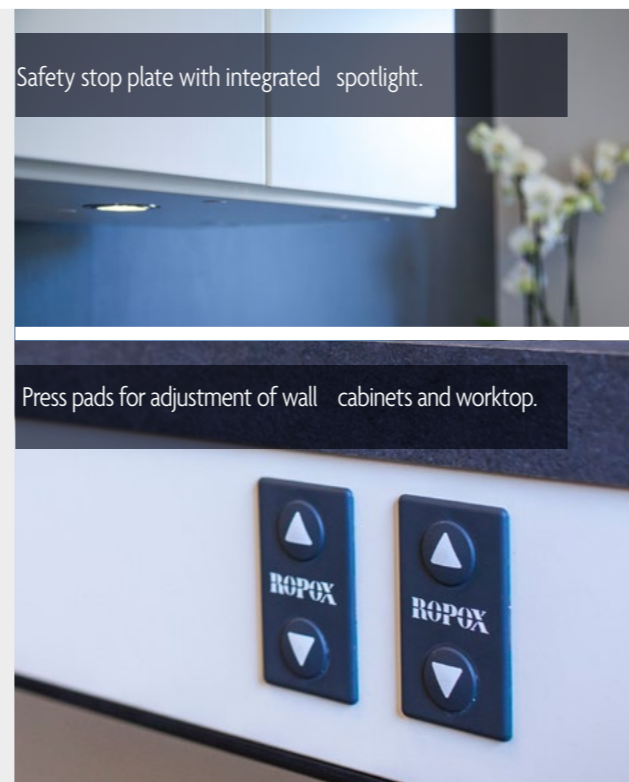

Safety stop plate with integrated spotlight.

Press pads for adjustment of wall cabinets and worktop.



## ROPOX Vertilinside (electric vertical height adjustment)

- Fitted inside the cabinet (standard widths between 400-1000mm).
- An excellent adaptation solution for existing units.
- Invisible from the outside.
- Adjustable down by 495mm.
- Designed to withstand loads up to 80kg.
- Height adjustable shelves.
- Lighting can be integrated in bottom safety stop plate.



## ROPOX Diagonal (electric diagonal height adjustment)


- Adjusts in a diagonal motion – 430mm down, 180mm forward.
- Frame depth only 115mm, hidden when cabinet is mounted.
- Frames available for cabinet widths from 400-1200,1300-1800mm.
- Designed to withstand loads up to 120/220kg.
- Lighting can be fitted in safety stop plate.



## ROPOX Brackets for wall cabinets


- Brackets enable stepped height adjustment of wall cabinets for very occasional adjustment. Cabinets are lifted off and placed at another height.
- Fitted on the back of the cabinets and therefore invisible.
- Height adjustment of cabinets in intervals of 42mm.
- Designed to withstand loads up to 90kg.

# HELP AND ADVICE ON KITCHEN LAYOUTS...

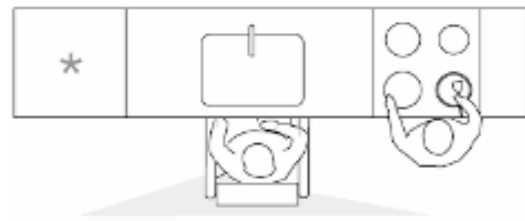

## ONE-SIDED KITCHEN

- Two people can work next to each other.
- Plenty of space for assistance.

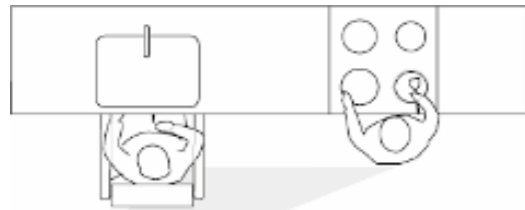


## TWO-SIDED KITCHEN

- People with mobility impairments can lean on worktops on two sides.
- Two independent work places.
- Sink and cooking section should be installed on same side.
- May be difficult for wheelchair turning and movement.




## L-SHAPED KITCHEN

- Seated and standing users could have independent work stations.
- Easy wheelchair access.
- Preparation area allows reach to sink and hob for easy movement of pans without moving position.

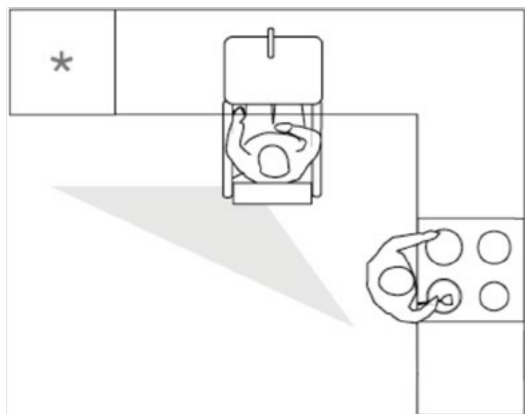


## U-SHAPED KITCHEN

- Objects can be pushed along the entire worktop.
- Easy wheelchair access.
- Greater distances to cover between appliances.

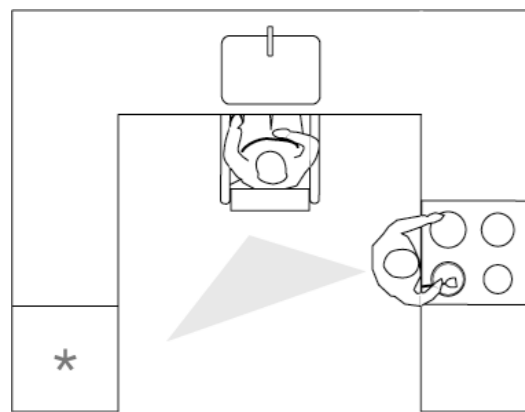


# OPTIMAL HEIGHT TO REACH RATIOS...

## ROPOX WALL CABINETS

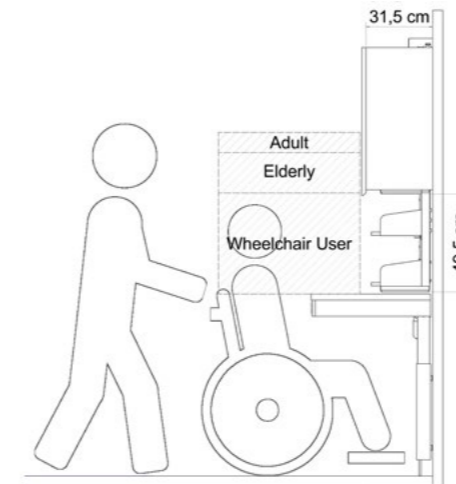

Ropox VertiInside

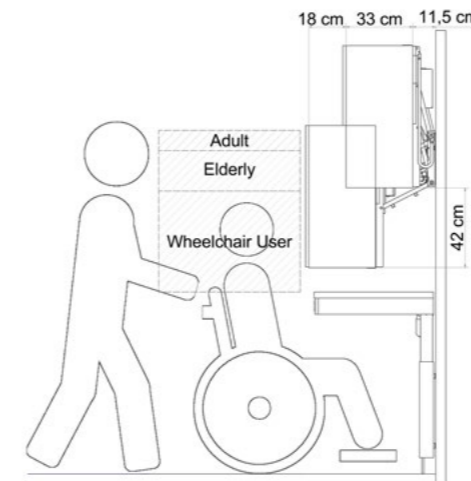

Ropox Diagonal

## OPTIMAL WORKTOP HEIGHTS

- Wheelchair users: 700-800mm
- Elderly: 850mm
- Adults: 900-950mm

# SAFETY IN FOCUS...



We recommend safety stop plates below all Ropox wall cabinets to avoid the risk of trapping. For further safety, wall cabinets can be connected to worktop frames via Ropox Smartboxes to avoid the trapping risk between wall cabinets and worktops, particularly when the worktops are moved up.



Ropox Smartboxes can avoid the trapping risk between height adjustable worktops installed next to each other. It's important that the side and front fascias have no worktop overhang.



Ropox safety strips ensure that worktops don't just stop, but reverse their movement immediately when activated.

# ABOUT ROPOX

## FREEDOM AND CONFIDENCE TO BE INDEPENDENT

“People with disabilities deserve an easier life with more self-esteem, quality of life and freedom to safely look after themselves.” This philosophy has been our goal since we set up our business back in 1962. Today, we are more than 60 employees who continue to pursue this mission.

With more than 60 years of experience, we are amongst the world’s leading specialists and producers of assistive devices and furniture for the people with disabilities as well as equipment for ergo- and physiotherapy and for rehabilitation. Being part of the AddLife group, we also have access to a large and strong network of partners who specialise in developing services and products to both the private and public healthcare sector.

We are represented in four continents of the world where our solutions are recognised and appreciated for their high quality, flexibility, and accessibility.

Find out more on [www.ropox.com](http://www.ropox.com).



### ROPOX TABLES

BEAUTIFUL HEIGHT-ADJUSTABLE TABLES FOR INDIVIDUAL NEEDS



### ROPOX KITCHENS

FLEXIBLE AND HEIGHT-ADJUSTABLE KITCHENS



### ROPOX BATHROOMS

INNOVATIVE AND EXCITING PRODUCTS IN A BEAUTIFUL COHESION

**ROPOX™**  
BETTER WAYS TO BETTER DAYS  
— A PART OF **AddLife**